

# ACTEON ENDODONTICS

**Endodontic Treatment Success  
with NEWTRON Devices and Tips**



# NEWTRON RANGE

**NEWTRON® DEVICES ARE VERSATILE,  
WITH THE MOST APPLICATIONS OF ANY ULTRASONIC SCALER**




## NEWTRON BOOSTER

Compact and efficient  
Non-LED Handpiece  
Quick-water disconnect



## NEWTRON P5 B0LED

LED Handpiece  
Quick-water disconnect

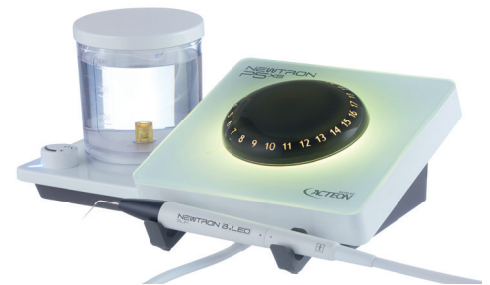


## NEWTRON P5XS B0LED


LED Handpiece  
Independent irrigation tank  
for any solution such as  
sodium hypochlorite or EDTA



PERIODONTICS

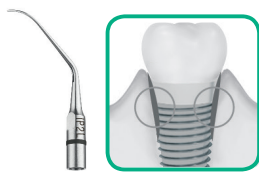

IMPLANT CARE

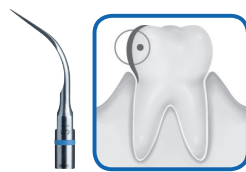

PROPHYLAXIS




ENDODONTICS



RESTORATION

BUILT  
IN +

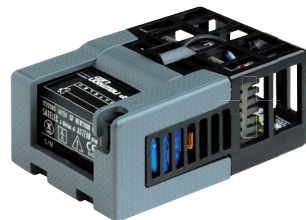


- Compact and lightweight module for easy installation
- Choice of several ergonomic handpieces that are easily disconnectable and autoclavable
- Cords are light and flexible, with perfect cord/handpiece balance and control during the procedure
- Extremely precise electronic power settings for unbeatable power and performance
- Can be integrated into any delivery system or dental chair

# ACTEON ENDODONTICS TIPS

## ENDOSUCCESS RETREATMENT KIT

Canal localization, access preparation, retreatment and loosening of posts or crowns


- ET18D and ETBD: cavity access finishing and calcified canal localization
- ET20: used in first coronal third extraction of filling material and removal of smear layer
- ET25 and ET25S - Titanium-Niobium: retreatment in the third root canal section and isthmuses. Corrosion resistant and can be pre-bent to treat curved canals, for removal of broken instruments and silver points
- ETPR: Post Removal

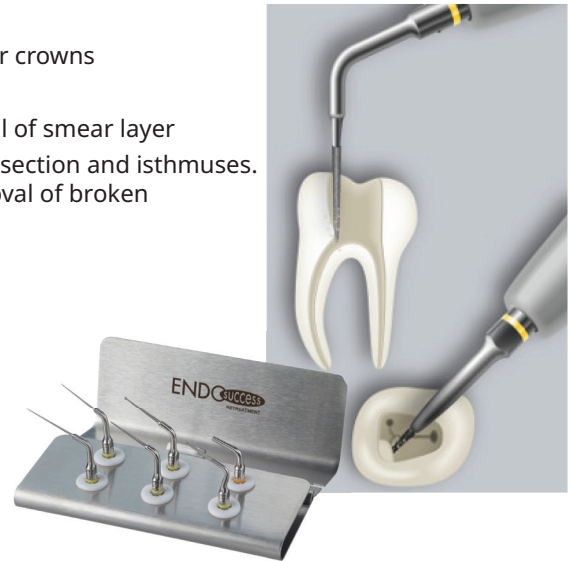

## CAP (CANAL ACCESS PREP) TIPS

Ultrasonic cavity access preparation: Perfectly-controlled ultrasonic vibrations for highly precise endodontic treatments

- 3 micro-blade tips that are the perfect solution for locating and opening hidden or calcified canals, shaping and finishing the access cavity
- The high strength of these tips decreases the risk of fractures, the rate at which they wear and increases the precision during procedures such as the shaping and widening of canal walls, the localization of canal orifices, and the fragmentation of calcifications and pulp stones
- CAP micro-blades tips as an alternative to diamond-coated tips

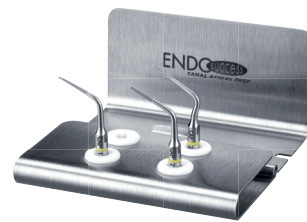

# ACTEON ENDODONTICS TIPS

## ENDOSUCCESS APICAL SURGERY KIT

5 ultrasonic instruments designed to address all anatomical configurations, preserving the suprastructures and root canals

- AS3D, AS6D & AS9D: 3-6-9 concept for easy and reproducible protocol in retrograde canal system preparation, up to 10mm
- Diamond-coated extremity (only): to avoid over-preparing the canal, for a more controlled retro endodontic treatment, to preserve bone and dental tissues



## IRRISAFE™

Passive ultrasonic irrigation: Completely safe and efficient smear layer, dental debris and bacteria removal in the root canal

- IrriSafe™ effectively cleans and disinfects root canals with hydrodynamic microcurrents
- IrriSafe™ guarantees a better prognosis for the non-surgical endodontic treatments

